

JASCO Spectra Manager – JASCO

V-730, V-750, V-760, V-770, V-780, MV-3000

INTRODUCTION

The Lumetics LINK™ software platform scans network locations for new measurement data files, copies data directly to a centralized database, and provides a powerful user interface for rapid multi-measurement multi-technique data aggregation, visualization, analysis, and reporting. LINK employs a client/server-based architecture where the LINK server hardware is provided by the end user and resides on the end user's network. The LINK client is a portable web-based application that may be placed on any computer with network connectivity to the LINK server. For successful import, the LINK webserver requires read access to the folders where user data resides.

The J-1000 Series Circular Dichroism spectrophotometers are specifically designed for high sensitivity measurements in the near- and far-UV regions. The J-1000 Series use high dynamic range detectors to obtain high quality spectra of samples with excellent signal-to-noise.

DETAILS

LINK requires either the far-UV or near-UV data set, exported from the JASCO Spectra Manager as TXT files. Optional support for multi-channel instrument.

The following raw data curves will be imported, in addition to all available instrument/analysis settings and parameters calculated by the instrument software:

- CD (mdeg) vs. Wavelength

Optional import data includes:

- Curve subtraction
- HT (V) as y-axis raw data, sub-populations based on HT ranges

Import Method options can be defined and calculated each measurement during import if it is first created/defined and then requested within an Import Method. The JASCO import method the option to specify an HT (V) acceptance range over which valid CD and HT raw data will be imported. Only data points corresponding to wavelengths where the HT (V) value is within this range will be imported. If these fields are left blank, LINK will import all CD and HT data found in the import file.

The LINKdb field, LNK_Import_HT(V)Range, will indicate if an import method has been applied.


EXAMPLES

Included below is a sample dashboard from JASCO measurement files:

1. Column Chart plotting raw data curves with visualization settings for "line + marker".



HT (V) & CD (mdeg) as a function of Wavelength

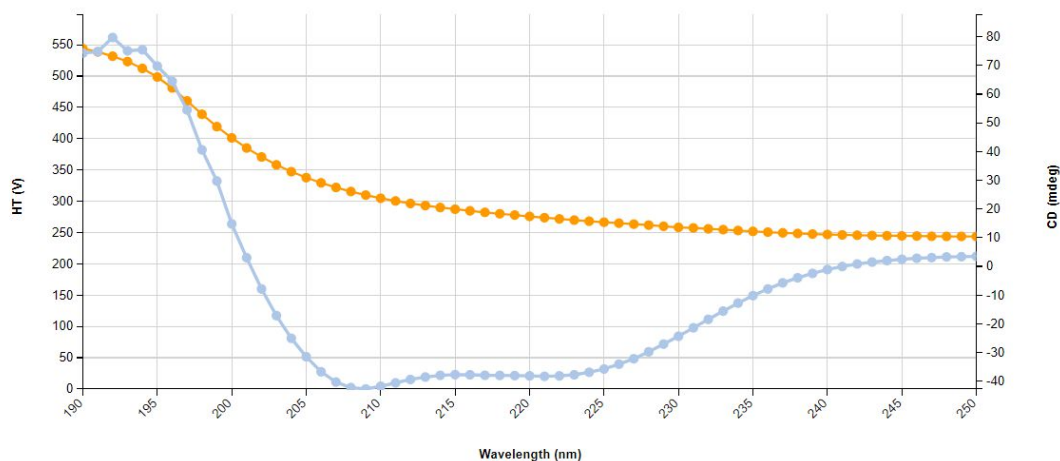

2. Tabular Summary examples

Measurement Summary Table – Measurement Results

| Sample Name | CD (mdeg) (Wavelength (nm) >1) - AVG | HT (V) (Wavelength (nm) >1) - AVG | Volume Analyzed (ml) - AVG | Temperature |
|-------------|---|--------------------------------------|----------------------------|-------------|
| Sample 1 | -1071.44 | 22758.84 | 10 | 25 |
| Sample 2 | -149.27 | 17701.28 | 15 | 30 |
| Sample 3 | 227.06 | 15677.38 | 20 | 18 |

Measurement Summary Table – Instrument Settings

| Sample Name | Scanning speed | Sensitivity | Data interval | Analysis Date | Time | Spectrometer/data system |
|-------------|----------------|-------------------------|---------------|---------------------|---------------------|--------------------------|
| Sample 1 | 10-8000 nm/min | -164 dBV re 1 µPa @ ... | 40 | 2018-09-28 16:23:47 | 16:23:47 | JASCO Corp., J-815, R... |
| Sample 2 | 10-8000 nm/min | -164 dBV re 1 µPa @ ... | 40 | Multiple (2 Values) | Multiple (2 Values) | JASCO Corp., J-815, R... |
| Sample 3 | 10-8000 nm/min | -164 dBV re 1 µPa @ ... | 40 | 2019-09-04 14:42:53 | 14:42:53 | JASCO Corp., J-815, R... |

DASHBOARD DOWNLOADS

Included below is a link to downloadable dashboards for JASCO measurement files:

<http://lumetics.com/dashboards/JASCO/JASCO.zip>

CONTACT LUMETICS

For direct assistance, please contact Lumetics LINK™ Support:

E-mail: support@lumetics.com

Phone: 1.613.614.874

Website: <http://lumetics.com/>

